


RealGUIDE™ Software Suite **Release 5.3**

OUT NOW

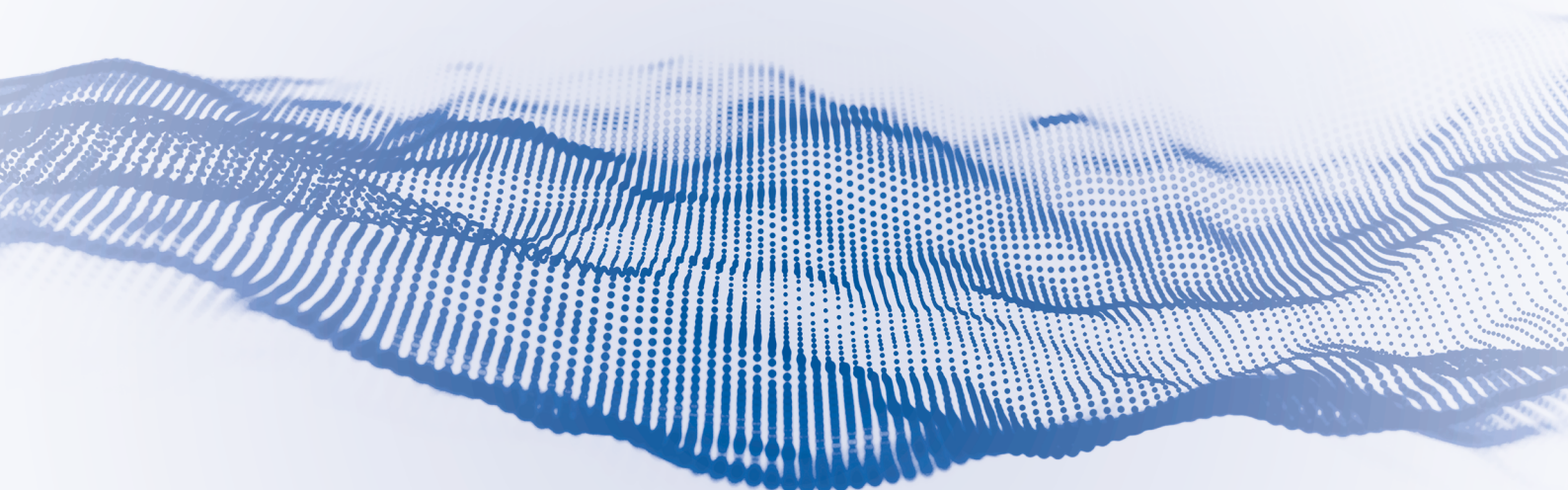


Evolving Dentistry Through Technology

RealGUIDE™ 5.3 Release Planner - New features


Assisted CPR drawing

Draw your patient's CPR with just a click!
Select the arch and RealGUIDE™ will draw the CPR for you.

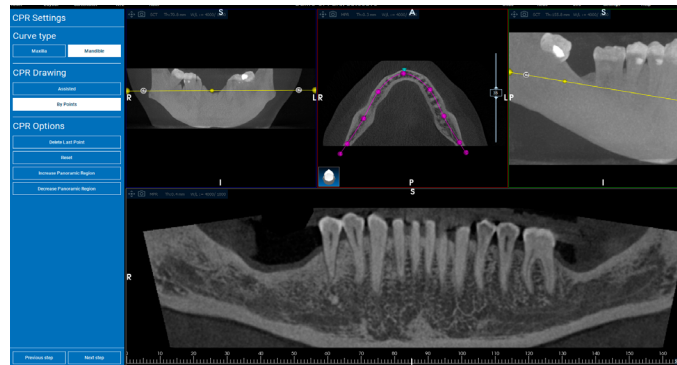


Advanced teeth setup

The new TEETH SETUP workflow is available in the PRO module and following. Design your wax-up from the teeth library in relation to the antagonist arch with a few clicks. You can also import your patient's photos for an aesthetic and functional design.

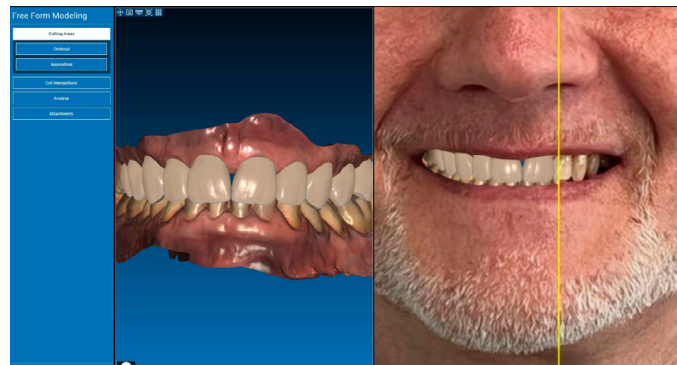


Assisted teeth segmentation

RealGUIDE™'s AI can help you with the most challenging operations! It can now segment the teeth visible in the CT exam and convert them to STL files.

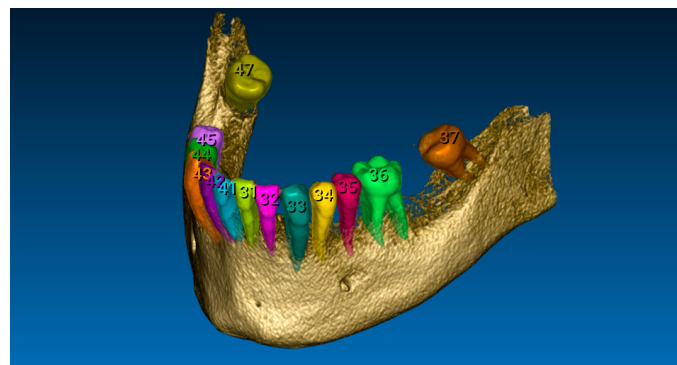

Assisted crowns segmentation

RealGUIDE's Artificial Intelligence is your best ally.
Perform partial or complete virtual extractions from the STL scan model effortlessly.

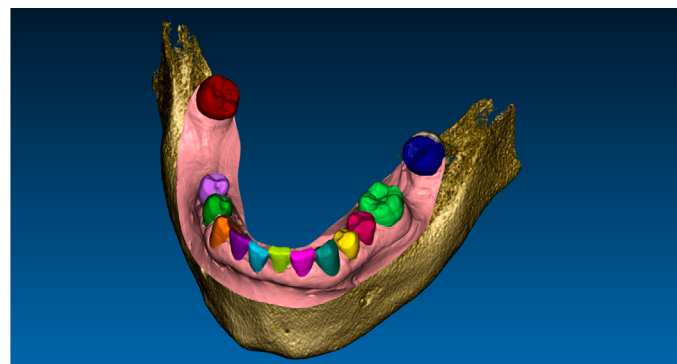


RealGUIDE™ 5.3 Release

CAD - New features


Occlusors library enhancement

Multiple item additions allowed for occlusors libraries

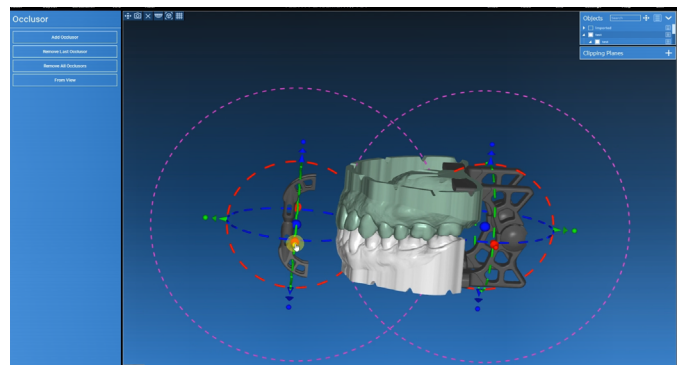


Smile Setup workflow

The new TEETH SETUP workflow is available in the PRO module and following. Design your wax-up from the teeth library in relation to the antagonist arch with a few clicks. You can also import your patient's photos for an aesthetic and functional design.

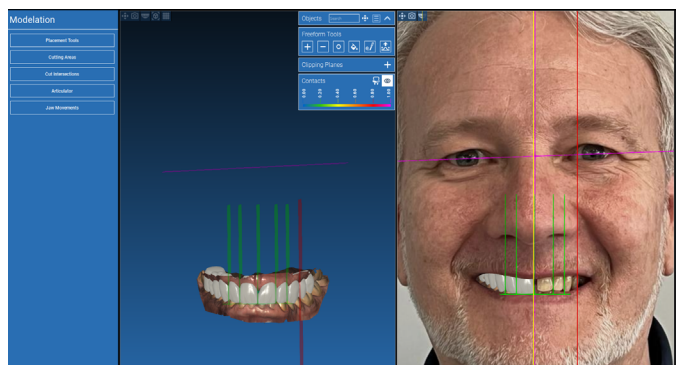


Faster immediate load workflows

Immediate loading on implants now even faster by machining simple cylinders.

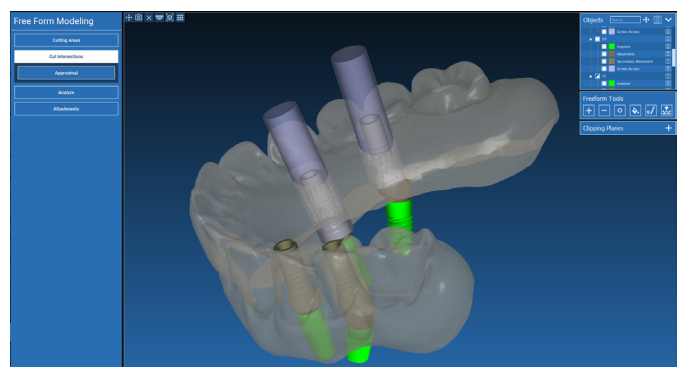


Addition of Ivoclar Stratos 300 virtual articulator

A new articulator is added to the many existing ones.

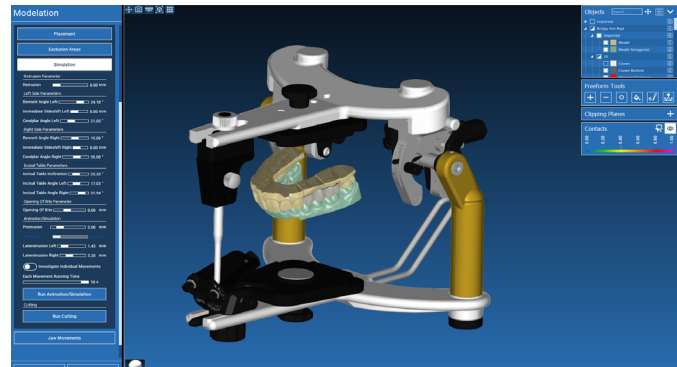


New files visualization adjustments

Adjusts the brightness and reflection of scans and modeling. Activates or deactivates texture greatly improving visibility during modeling.

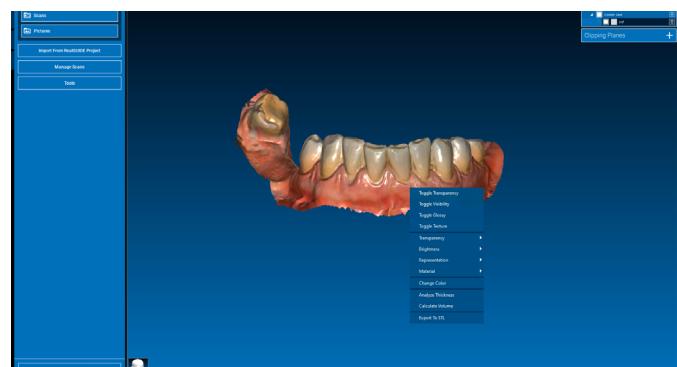

Other features:

- Smart Optics scanner integration
- Emergence profile control points addition/removing function
- Improved HyperDENT construction info export for complex reconstructions
- Improved HyperDENT construction info export for implant bars
- Reconstruction selection window further GUI optimization
- Teeth placement optimizations

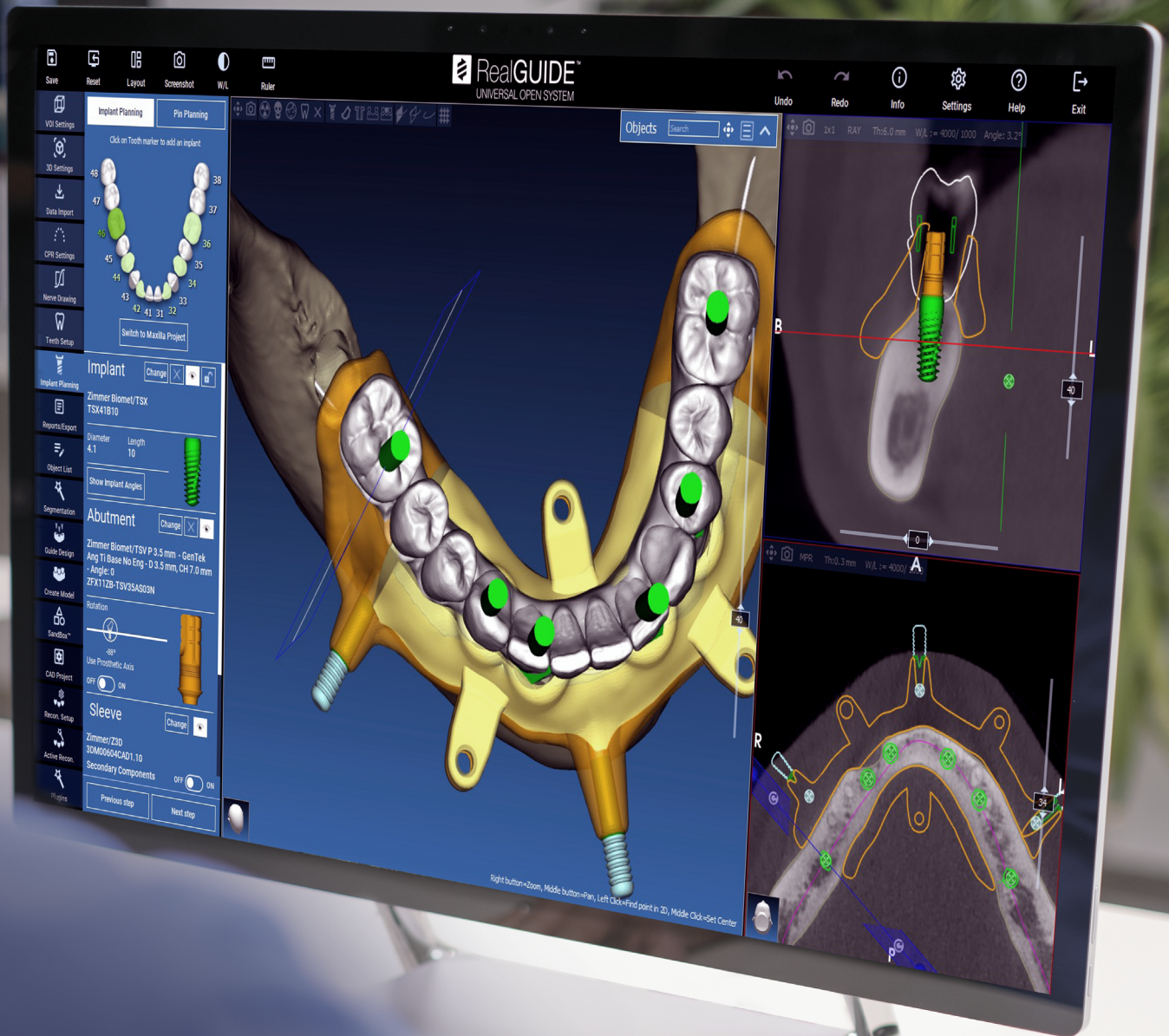

3DIEMME S.r.l.
Via Delle Baraggette 10 | 22060 Figino Serenza (CO) | Italia
Telefono +39 031.7073353
support@3diemme.it
www.3diemme.it

